


Intoximeters

Experience • Service • Integrity

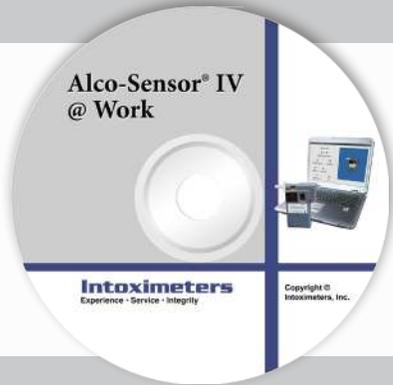

Alco-Sensor® IV @Work Software

Software application to manage an alcohol breath test program

Alco-Sensor IV @Work is a Windows® software application that combines the precision and accuracy of the Alco-Sensor IV with the processing power of your computer. This software package offers a low cost, easy to operate system for managing an alcohol breath test program (Alco-Sensor IV sold separately). Alco-Sensor IV @Work is designed to greatly simplify adhering to DOT testing protocol, in compliance with the Omnibus Transportation Employee Testing Act of 1991.

The @Work Software is an interactive interface which instructs the operator through a test sequence. Instructions for the sample collection process (DOT or non-DOT) are provided ensuring that the test is performed according to the programmed testing procedures and the Quality Assurance guidelines. The collected test information can be printed and is stored to the computer's hard drive for retrieval and report generation at a later time.

Intoximeters

Data Storage

Alco-Sensor IV @Work stores test results to a queryable database from which results may be reprinted and reports generated.

Data Collection

Alco-Sensor IV @Work eases and conduces uniform data entry allowing standard input data (i.e., technician or company information) to be stored to and picked from tables.

Printout Capability

Compatible, DOT approved (or look-alike) testing forms are used for test results.

Minimum Configuration Required to Run the Alco-Sensor IV @Work System:

Breath Analyzer

Intoximeters Alco-Sensor IV

PC Operating System

Microsoft Windows 7® or later.

Display (Monitor)

1024 x 768 (or better) True Color resolution

CPU, Available Hard Disk Space

Requirements do not exceed the minimum system requirements of the underlying Microsoft Windows® operating system.

Serial or USB Port

One available USB port or 9-pin serial port is required for connecting your PC to the Alco-Sensor IV. When using a USB port, a supported USB to serial adapter is required (included with package).

Communications (Cable)

Alco-Sensor IV @Work comes equipped with a 9-pin serial interface cable for the Alco-Sensor IV.

Parallel or USB Port

One available USB port or parallel port is required for connecting your PC to a supported printer.

Test Result Printer*

A supported parallel or USB printer (purchased separately). Choose either a supported label printer (you affix labels to appropriate sections of each required copy of the testing form) or a supported dot matrix impact printer (prints directly onto the carbonless triplicate testing form).


Report Printer (optional)

Summary reports and statistics may be printed to various laser or inkjet printers supported by Windows®.

* Please visit <http://www.intox.com/@Work> for the most current list of supported full form and label printers.

Alco-Sensor® IV @Work Software

Software application to manage an alcohol breath testing program

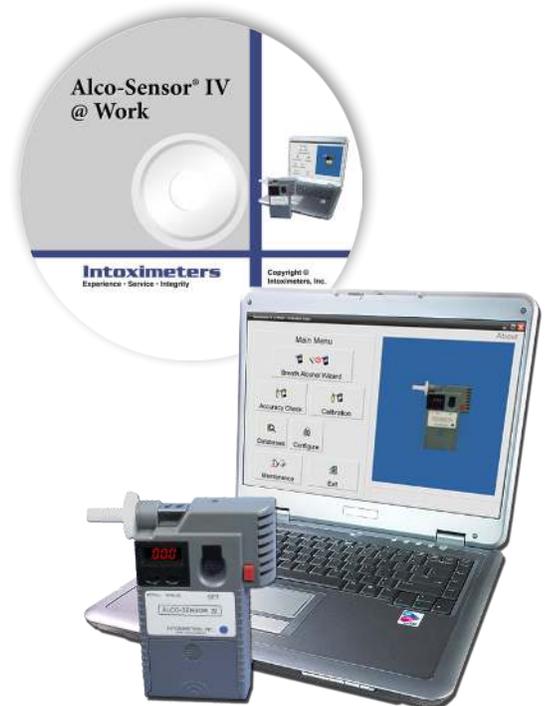

Microsoft® is a registered trademark of Microsoft® Corp. in the United States and/or other countries.

The names of actual companies and products mentioned herein may be the trademarks of their respective owners.